
Remote Microphone Connection Procedures for students with Oticon rechargeable hearing aids


– Roger Touchscreen as Transmitter & EduMic as Intermediary

Initial Set up:


1. Pair the EduMic to the hearing aids (to be done during initial set up only).
   a. Turn OFF the hearing aids by pushing and holding the LOWER button for 3 seconds until indicator light flashes orange.

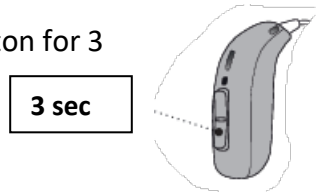

   b. Turn ON the EduMic by pressing the POWER button for 2 seconds until both indicators turn WHITE then release the button. Wait 10 seconds until the indicator turns GREEN.

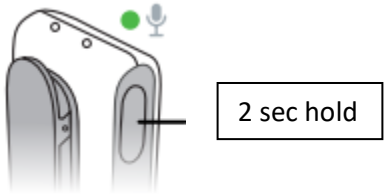

   c. Power ON the hearing aids by pushing and holding the LOWER button for 2 seconds until light flashes green.

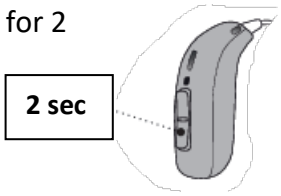

   d. Hold the hearing aids within **8 inches** of the EduMic. EduMic flashes BLUE while pairing is in process. This may take up to 30 seconds. When pairing is complete, light turns GREEN.

2. **CONNECT** the Roger X receiver to the Roger Touchscreen transmitter. (**CONNECT** needs to be done at the initial set up only; connection holds until the Roger X is connected to another Roger transmitter.)

   a. Turn ON the Roger Touchscreen by pressing and holding the side power button for 3 seconds.



b. Insert the Roger X into the FM connector of the powered-on EduMic.



c. While holding the powered-up Roger X within 1 inch of the Roger Touchscreen, tap CONNECT.



d. If the operation is successful, the message “Roger X is connected” will appear on the Roger Touchscreen.

<u>Daily Set Up</u>:

At the beginning of the school day:

a) Wait until the Hearing Aid user is in the same room as the Teacher. Teacher wears the Touchscreen using the supplied lanyard so that the top of Touchscreen is within 20 cm of mouth.



b) Student put the EduMic around their neck using the supplied lanyard.

c) Power ON the EduMic (with Roger X inserted).

d) Conduct functional check by giving the student a simple, unexpected, instruction from a distance of at least 3 metres with no visual cues (turn your back or cover the lower face). E.g. “Can you point to the door?” or “Touch your head”. Verify student’s understanding.

During Instruction Breaks (e.g. individual work time, recess):

a) Teacher to mute the Roger Touchscreen by pressing the Mute button once. A RED light will blink around the MUTE button.




b) When instruction resumes, unmute the Roger Touchscreen by pressing the Mute button once.

At the end of the school day:

a) Turn OFF the Roger Touchscreen by pressing and holding the power button for 3 sec and tap on the upper right check mark for Yes.




b) Turn OFF the EduMic by pressing and holding the POWER button for 2 seconds until the indicator light turns off.

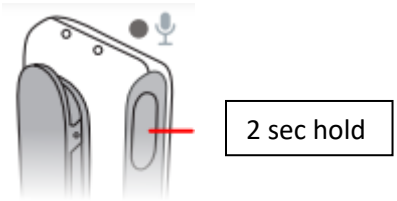

c) Charge both Touchscreen and EduMic in a safe, dry, and clean place.



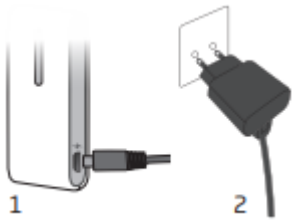

<u>For storage before a weekend or school breaks</u>:

- ➔ Detach the Roger X receiver from the EduMic to allow receiver contacts to air dry.
- ➔ Store Roger X receiver in a secure, dry place.

For equipment support, contact Auditory Outreach audiology at 1-866-430-4327, or email carrie.siu@sd47.bc.ca .