

Easy Sheet for a Roger Neckloop + Oticon Hearing Aid set up for Speech-to-Text

Roger DM:

Pair the EduMic to the Hearing Aids (HA). (This needs to be done at the initial set up only. Pairing holds until EduMic is manually unpaired from HA). Refer to easy sheet for Oticon Hearing Aid with EduMic for steps on how to pair. Leave edumic turned on.

Connect Roger X receiver to bottom of Edumic.

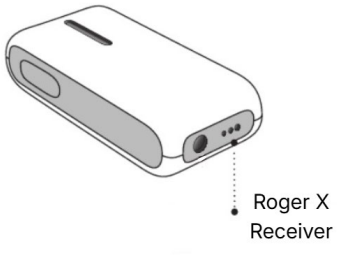

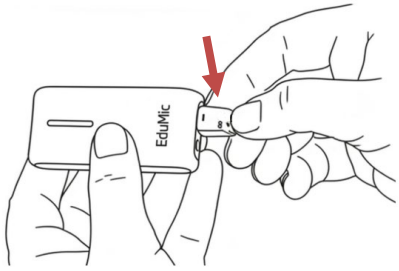


Connect Roger Touchscreen or Roger ON to Roger X receiver. Need to be within *3 inches* of the receiver to connect.

Roger Touchscreen:



Roger ON:

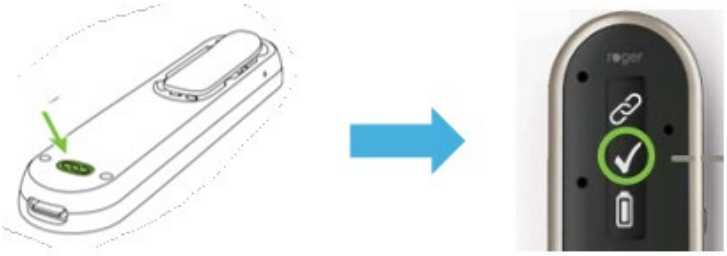

Student's Oticon hearing aids are now connected to a Roger DM.

Speech-to-Text:

Connect Roger Device to Roger Neckloop.


Initial set up:

1) **Power ON** the Neckloop by sliding the toggle up.



2) **Connect** the Neckloop to the transmitter by holding transmitter within 4 inches of the Neckloop and tapping "Connect". For a Roger ON, press the connect button on the back.

Touchscreen:

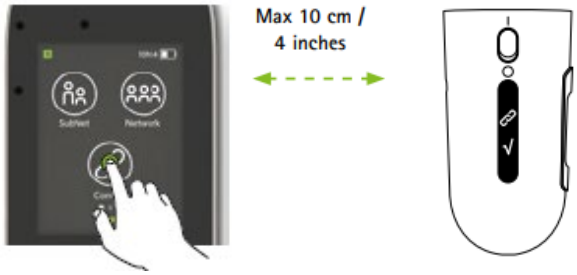


Roger ON:



Note: For a Roger touchscreen, if connection is successful, *"Roger Rx is connected"* will appear briefly on the Touchscreen. For a Roger ON, if connection is successful, you will see a checkmark on the screen.

3) Using the **USB-C** cord and adaptors as necessary, connect Roger Neckloop to student's personal device (laptop, phone, or iPad) and open speech-to-text software. (Please refer to Phonak STT Roger Neckloop Guide on AO website for info on adaptor cables).

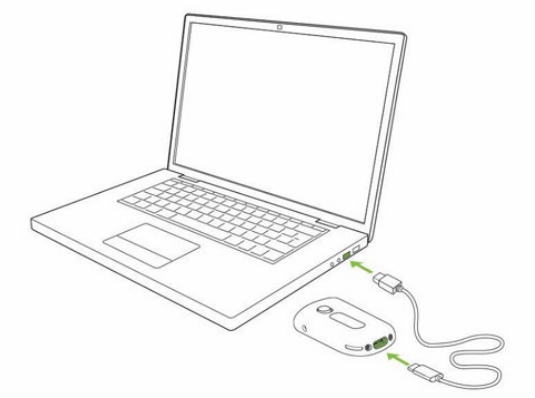

Teacher can now wear the Roger Touchscreen/ON. Student wears edumic or places it on their desk.

Charge all devices at the end of each day.

For tech support, contact Auditory Outreach at 1-866-430-4327, info@auditoryoutreach.ca